





Influenza H7 Pseudovirus

Catalog #: 1005-PP-007

Description:

The influenza H7 pseudovirus is a third-generation lentiviral vector pseudotyped with the hemagglutinin protein from the influenza A/Shanghai/4664T/2013 (H7N9) strain. It contains a firefly luciferase reporter system. This pseudovirus is classified under biosafety level 2 (BSL-2) and is not known or suspected to contain any replicative biological agents due to multiple safety modifications in the viral genome (e.g.: altered 3' LTR rendering the vector self-inactivating, exclusion of the viral tat gene, and replacement of the U3 region in the 5' LTR with a strong tat-independent constitutive promoter).

Formulation & Storage: Supplied in DMEM supplemented with 10% Fetal Bovine Serum. Store at -80°C. Avoid Thaw/freeze cycle.

Notes & Usage Guidelines:

Recommended Dilution: 1/200 to 1/100 **Luciferase Units for Assays:** 100,000 to 220,000


Reporter: Luciferase

~No. of plates that can be tested with a single aliquot: 1 plate

Applications

Viral Titration:

Pseudoviruses were serially diluted 2-fold in DMEM supplemented with 2% FBS, 25 mM HEPES, and 1x GlutaMAX (transduction medium). Next, 50 μ L of transduction medium containing 10,000 HEK293-T cells was added to each well. Twenty-four hours post-transduction, 100 μ L of DMEM supplemented with 10% FBS, 25 mM HEPES, and 1x GlutaMAX was added. Luciferase expression was measured 72 hours post-transduction using the Bright-GloTM Luciferase reagent (Promega, Madison, WI, USA) according to the manufacturer's instructions. Analysis was performed using a Biotek Synergy H1 plate reader (Gain: 135). Cells with media alone were used as a control. Data are representative of duplicates. Viral titer was expressed as Relative Luciferase Units per microliter (RLU/ μ L). The titration obtained may vary depending on the user's instrument, handling, or luciferase assay system.

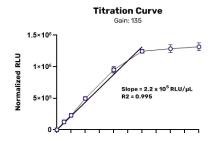

Certificate of analysis:

A hardcopy of datasheet is sent along with the products. Please refer to it for detailed information. For older lots, refer to the applicable certificate of analysis that may be requested at services@ibtbioservices.com

Related Products:

IBT provides a wide array of anti-filovirus specific antibodies, recombinant proteins, and other infectious disease reagents. Please see our website, www.ibtbioservices.com for more details.